

TISSHIELD™

Injuries can be avoided




SHIELD
Your Patient
Your Team
Yourself




Retracting & Shielding Liver
from Thermal Damage
(Liver Resection)

**FROM
EXPOSURE
TO
CLOSURE**

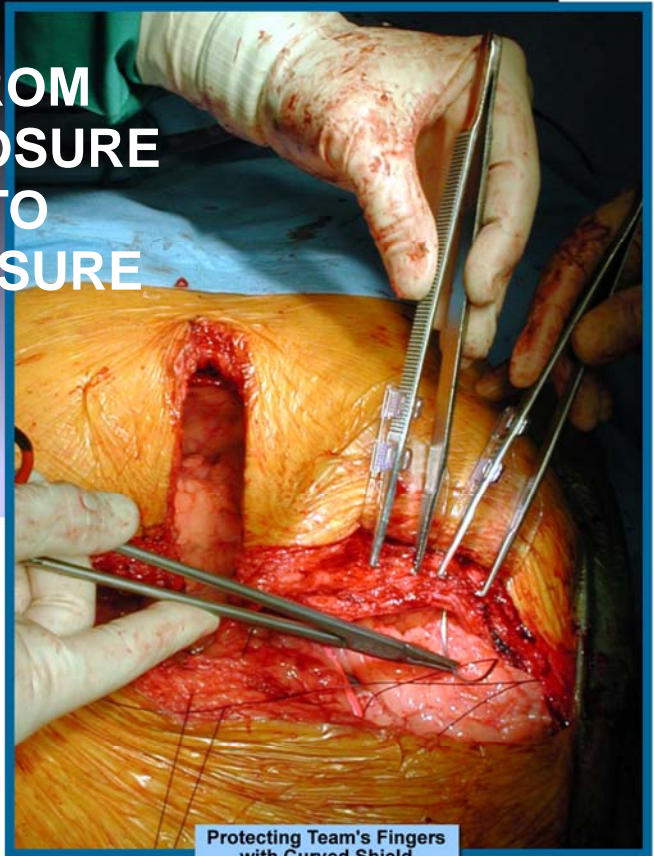

Protecting Team's Fingers
with Curved Shield
(Approximating Fascia)

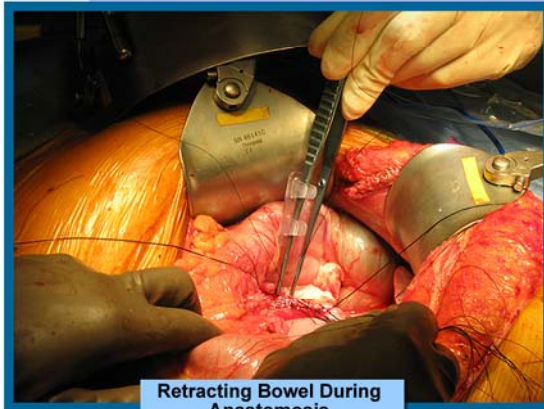

Retracting Bowel During
Anastomosis
(Gastrojejunostomy)

INTRODUCING TISSHIELD

APPLICATIONS:

Injuries can be avoided. The Tisshield™ is a versatile forceps attachment that helps surgeons and their operating teams avoid percutaneous injury, increases surgical efficiency, and protects patients' tissues. Tisshield™ prevents injuries by keeping fingers out of harm's way, while shielding patients' tissues against inadvertent mechanical or thermal damage. The device also allows forceps to serve as transparent retractors, gently providing optimal site exposure in tight quarters while retaining obstructive surrounding tissues. Above all, the Tisshield™ is inexpensive, easy to use, and is recommended by leading surgeons.

*Shield your patient. Shield your team. Shield yourself.
Use Tisshield™, the Essential Forceps Companion.*

- Procedures where exposure is essential.
- Procedures involving incisions greater than the device width (1.25 in. or 3.175 cm).
- Surgical sites involving tissue or organs large enough to be accommodated by the device, but not requiring more than 3lbs of force to displace such tissue.

USAGE:

- Use the Tisshield™ to retract skin, tissues, and organs during suturing or handling of tissues. Avoid using fingers whenever possible.
- Use Tisshield™ to shield the patient's tissues against inadvertent mechanical or thermal damage.
- For increased safety, surgical team members can each use Tisshield™ to perform retraction.

FEATURES AND BENEFITS	FORCEPS	FINGERS	TISSHIELD™
Increases view of surgical site?	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Enables suturing without an assistant?			<input checked="" type="checkbox"/>
Avoids finger exposure?			<input checked="" type="checkbox"/>
Prevents tissues from impeding surgical site?		<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Retracted tissues are clearly visible?	<input checked="" type="checkbox"/>		<input checked="" type="checkbox"/>
Uses wide surface area to minimize pressure against tissues?		<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Shields tissues against mechanical damage?			<input checked="" type="checkbox"/>
Protects tissues against thermal damage?			<input checked="" type="checkbox"/>

WHAT LEADING SURGEONS HAVE TO SAY ABOUT TISSHIELD™:

Tisshield™ is effective as a surgical tool, overall. I would recommend Tisshield™ to others.

-Sean C., M.D.
Hepatobiliary & Pancreatic Surgery
(America's Best Doctors, 2004)

Tisshield™ enhances the effectiveness of forceps. Tisshield™ improves the speed of surgical operations.

-Tam L., M.D.
Colon & Rectal Surgery
(2004 Physicians of Excellence)

ORDER NUMBER	ITEM DESCRIPTION	SPECIFICATIONS
TS1002	Box of 10 Sterilized Tisshield™ Kit P/N TS1001. Two (2) TISSHIELDS™ are provided in each kit P/N TS1001: The straight shield P/N TS1100 and the curved shield P/N TS1200	APPLICABLE FORCEPS: Tissue / DeBakey FORCEPS LENGTHS: 15cm-24cm FORCEPS JAW SIZES: 1,5mm-9mm Width, 2,2mm-3mm Thick MATERIAL: GE Advanced Thermoplastic FDA/USP Class VI Compliant. BLADE SIZE: 3,2cm Width X 5,7cm Length: STERILITY: Gamma Radiation Sterilized WEIGHT: 4,5 grams ea. USAGE: SINGLE USE ONLY



TO PLACE AN ORDER, CALL 1-877-TISSHIELD (1-877-847-7443), OR VISIT www.tisshield.com

Dexteus
Innovations for Better Surgery

P.O. Box 51810
Irvine, CA 92619
USA
Phone: 1-866-4DEXTEUS (866-433-9838)
Fax: 949-305-9784